
TESLA HEATING MASSAGER


USER MANUAL

- EN
- DE
- CZ
- SK
- PL
- HU
- RO

DEVICE DESCRIPTION

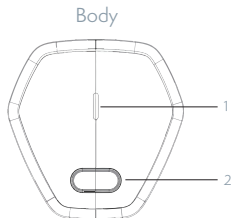- 1) Status indicator (working/charging)
- 2) Power button

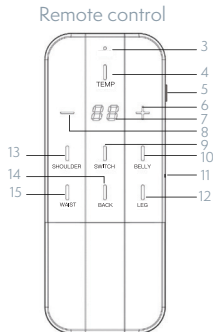- 3) Battery/charging indicator
- 4) Temperature setting button (green/yellow/red)
- 5) Power button / Lock button
- 6) Intensity increase button
- 7) Digital display for intensity
- 8) Intensity decrease button
- 9) Switch button
(Click to switch between parts of the body)
- 10) Abdomen massage indicator
- 11) Reset button
- 12) Leg massage indicator
- 13) Shoulder massage indicator
- 14) Back massage indicator
- 15) Waist massage indicator

INSTRUCTIONS

1. Installation

A) Remove the blue gel protective film and attach the gel pads to the thin electrical coating on the back of the device.

- Before use, wipe the electrical coating carefully to prevent poor adhesion or disturbance of electrical conductivity caused by body sweat and grease.

- Make sure that the gel pads and the electrical coating fit tightly without any bubbles. If bubbles appear, use your finger to press the gel pads, which are separated by the transparent film.

B) Install the body with correctly determined left and right direction of wing extension. (Make sure the magnet is aligned).

2. Power on

Powering on for the first time (pairing the device with remote control):

Long press the power button on the remote control and the body for 6 seconds while the body and the remote are both off. The body's white light and the remote's switch light will begin to flash. When the white light on the body remains on and the remote switch light changes to white and remains on, pairing is successful.

Normally powering on:

Long press the power button for 3 seconds while the body and the remote control are both off. When the device is successfully connected to the remote control, the body's white light, digital display, body area indicator, and intensity decrease and increase buttons will remain on.

Pairing is required the first time the product is used. The body will automatically pair with the remote control after normal power-up.

3. Attaching the gel

Remove the white protective film of the gel and stick the product on the place where the massage is desired.

4. Start massage

Press the „+“ button on the remote control. The digital display screen will change from 00 to 01 and the body will start the massage.

5.Intensity

Press the „+/-“ button to adjust the intensity and select the intensity most suitable for you. The intensity is adjustable from 0 to 15.

- Long press the „-“ button to switch the massage output.

6.Mode

Press the switch button on the remote control to switch between massaging different parts of the body. There are 5 massage modes for different parts of the body.

7.Heating function

Press the temperature setting button to set the heating function. The heating temperature is divided into three levels: low, medium, and high.

- After 8 seconds of sticking it to the body, reheat it for a better experience.

8. Child lock function for remote control

When not operated, the remote control automatically enters the child lock function after 30 s. When you need to use the remote control, you can wake it up by pressing the switch button, and you can also press the switch button to lock the remote control screen during use.

9. Shutdown

The device switches off automatically after 15 minutes or by long pressing the power button on the remote control for 3 minutes.

10. Storage

After use, conveniently stick the white protective film back onto the gel and store it in a cool, dry place.

MAINTENANCE

1. To prolong the life of the built-in battery, store the massager when fully charged when not in use for an extended period of time and charge it once every 2-3 months.
2. If dust or sweat appears on the surface of the massager, wipe it gently with damp wipes (the wipes should not leak water).
3. If the gel pads are difficult to attach due to sweat, follow the steps below:
 - Use your finger to scoop up a drop of water and apply it evenly to the gel or gently wipe away dust.
 - After thoroughly wiping off the water, place the gel on top and place it in a well-ventilated and shaded area to dry.
 - Air dry until no moisture remains on the surface, stick the gel pads onto the white protective film, store them in a cool, dry place. (If the stickiness of the gel cannot be restored, please replace it with a new one.)

4. Method of gel removal

- Pinch one end of the gel pads with your fingers and tear off. (It is better to roll them up when tearing them off.)
- Completely remove the residual gel pads on the conductive coating to avoid the negative impact of incomplete cleaning. (Can be more easily removed by sticking a small amount of residual gel with a torn gel pad.)
- Make sure the machine is switched off when removing the gel pads, otherwise you may get burnt or scalded.
- After peeling off the gel, wipe the conductive coating carefully to prevent lack of adhesion or conductivity due to greasy sweat or body perspiration or grease.

5. Do not use or store this product in the following places:

- Direct sunlight
- Near fireworks
- Hot and humid place
- Dusty area
- Place with water
- Place with strong electromagnetic radiation

Suggested massage parts

Shoulder

Back

Abdomen

Waist

Leg

Not suggested massage parts

Head

Neck

Chest

Sensitive areas

TROUBLESHOOTING

Problem	Possible causes	Solution
The device will not switch on when the ON/OFF button is pressed.	The battery may be low.	Fully charge the battery. If the device still cannot be switched on, contact customer service.
The device does not stick to the body.	Gel pads may be dirty or worn; the skin area to be treated may not be clean enough.	<ol style="list-style-type: none">1. Clean the adhesive surface of the gel pads with a damp, lint-free cloth. If the gel pads still do not hold firmly, they need to be replaced.2. Cleanse the skin before any application; do not use skin care creams or oils before treatment.
There is no noticeable stimulation.	Gel pads may be dirty or worn. Batteries may need charging.	<ol style="list-style-type: none">1. Make sure the gel pads are in firm contact with the area to be treated.2. The battery is nearly empty. Fully charge it.
An unpleasant sensation appears at the site of application of gel pads.	The gel pads are worn. This can cause skin irritation.	Replace gel pads.

<p>The skin in the treated area turns red.</p>	<p>If the redness under the gel pads disappears quickly, this is not dangerous and is due to locally stimulated and increased blood flow.</p>	<p>Stop the process immediately and wait until your skin returns to normal. However, if skin irritation persists and is accompanied by itching or inflammation, consult your doctor before continuing treatment as this may be due to an allergic reaction to the sticky surface.</p>
<p>The skin turns red under the heat.</p>	<p>Redness under the thermal surface of the massage paste is due to the strongly stimulated blood flow in the heated tissue.</p>	<p>The skin should return to normal a few minutes after treatment. If redness persists and if itching or inflammation occurs, consult your doctor before using the heat function again.</p>
<p>The massage paste is too hot.</p>		<p>Press the ON/OFF button to stop the heating function and continue massaging.</p>

WARNINGS AND SAFETY INSTRUCTIONS

Contraindications

1. People with implanted electrical devices.
2. People with known or acute cardiac arrhythmia or cardiac impulse and conduction system disorders.
3. Pregnant women and children under 12 years of age.
4. All types of pets.

Persons with medical conditions on the advice of a doctor

1. Persons who are being treated or taking medication for a serious illness.
2. Persons with fever.
3. Persons with acute inflammation, purulent inflammation, and tendency to bleed.
4. Persons with blood abnormalities, infectious diseases, or skin diseases.
5. Persons with skin allergies.

Notes on use

1. Keep out of reach of children.
2. Elderly persons and persons who cannot take care of themselves should be accompanied during use.
3. Do not put the product or gel pads in the mouth.
4. Do not use this product during strenuous exercise.
5. Keep skin dry during use.
6. Do not apply to a wound.
7. When using, please avoid metal objects such as necklaces, bracelets, etc.
8. If you feel uncomfortable during use, please stop using.
9. Do not use continuously for too long (recommended no longer than 30 minutes) to prevent muscle fatigue, skin tingling, or redness.

Technical Specifications

Material: Silicone + ABS + PC

Size: 274 x 75 x 22 mm

Weight: 140 g

Rated voltage: 5 V

Rated power: 7.5 W

Battery capacity: 500 mAh

Charging voltage: DC 5 V / 3 A

Battery charging time: about 2 h

(This product is not supported by an external mobile power supply.)

Package contents:

1x Body

1x Remote control

1x Left wing

1x Right wing

1x USB charging cable

4x Gel pads

User Manual

The package of this product does not include a power adapter. The power adapters applicable to this product are power adapters for various mobile phones and computer USB interfaces. Using other types of power adapters may damage the product.

Light indicator description

Body lights status:

- White light is on: the device is in the on/standby mode.
- White light is flashing: the device is pairing with the remote control.
- Indicator light gradient: the unit is operating independently.
- One of the colours (red/yellow/green corresponds to low/medium/high temperature) remains lit: the device is warming up independently.
- One of the red/yellow/green gradient colours remains lit: the unit and the heating mode are working together.
- Blue light flashes: low power/charging status.

Status of the remote control indicator light

- The digital display and the body part light stay on: the device is powered on.
- Switch button is flashing: the device is in the Bluetooth pairing state.
- Digital display >00: Heating light does not light and body part light remains lit: TENS state is operating independently.
- One of the colors red/yellow/green (corresponding to low, medium, high temperature) remains lit: the device warms up on its own.
- One of the colours (red/yellow/green corresponds to low/medium/high temperature) remains lit: EMS and heating are working together.
- Red light flashing: the device needs charging.
- Blue light flashing: the device is charging.
- Blue light remains lit: the device is charged.

INFORMATION ABOUT DISPOSAL AND RECYCLING

This product is marked with the symbol for separate collection. The product must be disposed of in accordance with the regulations for disposal of electrical and electronic equipment (Directive 2012/19/EU on waste electrical and electronic equipment). Disposal together with regular municipal waste is prohibited. Dispose of all electrical and electronic products in accordance with all local and European regulations at the designated collection points which hold the appropriate authorisation and certification in line with the local and legislative regulations. Correct disposal and recycling help to minimise impacts on the environment and human health. Further information regarding disposal can be obtained from the vendor, authorised service centre or local authorities.

EU DECLARATION OF CONFORMITY

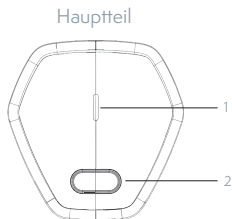
Hereby, Tesla Solar s.r.o declares that the radio equipment type TSL-HC-HM is in compliance with Directive 2014/53/EU. The full text of the EU declaration of conformity is available at the following internet address: tsl.sh/doc

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